

Medtronic

MiniMed™

**Quick-set™ Infusion Set and Reservoir, and Guardian™ 4
Sensor**

Quick Reference Guide



MiniMed and MiniMed logo are trademarks of Medtronic MiniMed, Inc. ^{TM*} Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

Contents

- **Changing the MiniMed™ Quick-set™ Infusion Set and Reservoir 7**
The MiniMed Quick-set infusion set and the reservoir can be used for up to three days. In most cases, the infusion set and reservoir should be changed at the same time.
- **Disconnecting and Reconnecting the MiniMed Quick-set Infusion Set 23**
- **Changing the MiniMed Reservoir Only 27**
The MiniMed reservoir can be used for up to three days. If the reservoir runs out of insulin and the infusion set has not been used for three full days, the New Reservoir Only option may be used to change the reservoir. The infusion set must still be replaced within three days of use.
- **Changing the MiniMed Quick-set Infusion Set Only 39**
The MiniMed Quick-set infusion set can be used for up to three days. If only the infusion set needs to be changed before three full days of use, the New Set Only option may be used to change the infusion set. The reservoir must still be replaced within three days of use.
- **Inserting the Guardian™ 4 Sensor with One-press Serter 55**

WARNING: The following warning applies to changing the reservoir and infusion set, changing the reservoir only, and changing the infusion set only.

Always confirm that the infusion set tubing is disconnected from the body before doing the following steps:

- placing the reservoir into the pump
- rewinding the pump
- loading the reservoir
- filling the infusion set tubing

Failing to disconnect the infusion set tubing from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

Changing the MiniMed™ Quick-set™ Infusion Set and Reservoir



The MiniMed Quick-set infusion set and the reservoir can be used for up to three days. In most cases, the infusion set and reservoir should be changed at the same time.

Note: The MiniMed Quick-set infusion set and reservoir can be used for up to three days. To avoid wasting insulin, fill the reservoir with enough insulin to last until the next infusion set change.

START HERE

Follow these steps to change the infusion set and reservoir.

1




Wash hands with soap and water. On the pump, press  to go to the Menu screen.

2



Select , and then select **New Reservoir & Set**.

3

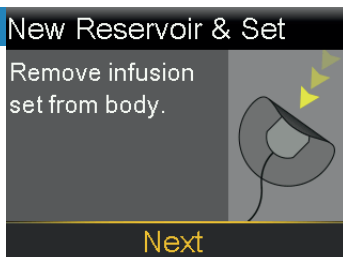


Remove the infusion set by loosening the adhesive and pulling the set away from the body. Select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before rewinding the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

4 New Reservoir & Set
 Remove reservoir from pump.
 Then select Rewind.

Rewind

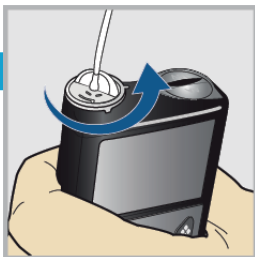

Remove the used reservoir from the pump.

5 New Reservoir & Set
 Remove reservoir from pump.
 Then select Rewind.

Rewind

Select **Rewind**.
Do not connect the infusion set to the body.




6 New Reservoir & Set
 1. Fill new reservoir.
 2. Connect new tubing to reservoir.
 DO NOT CONNECT TO BODY.


Next

Follow the next steps to fill the new reservoir with insulin and to connect the new infusion set tubing.
 Do not select **Next**.


FILL THE RESERVOIR

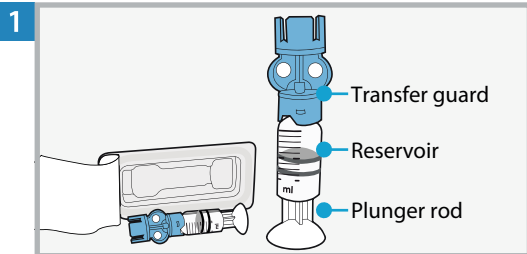


Remove the reservoir from the package. Make sure the insulin vial is at room temperature to reduce the risk of air bubbles.

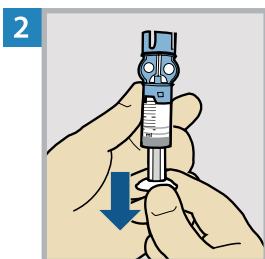

Pull the plunger down based on the planned 2-3 day insulin fill amount.


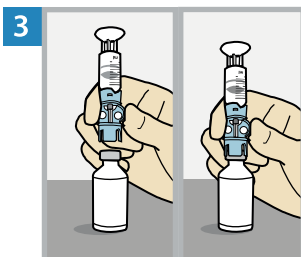


Wipe the top of the vial with alcohol. Place the vial on a sturdy, flat surface. Firmly press the transfer guard onto the vial.

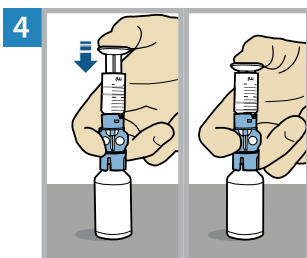


Push and hold the plunger down.


WARNING: Do not use the reservoir or infusion set if insulin or any liquid gets on the top of the reservoir or inside the tubing connector, as shown in the image. Insulin or any liquid can temporarily block the vents. This may result in the delivery of too little or too much insulin, which can cause hyperglycemia or hypoglycemia. If insulin or any liquid gets on the top of the reservoir or inside the tubing connector, start over with a new reservoir and infusion set.

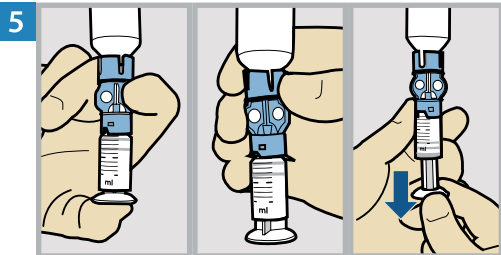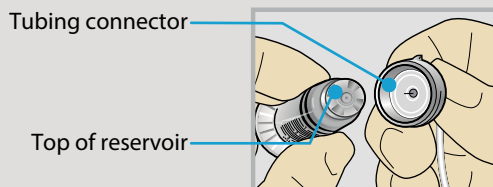

Keeping a thumb on the plunger, flip the vial over so the vial is on top. Release the thumb and pull the plunger down to fill the reservoir with insulin.

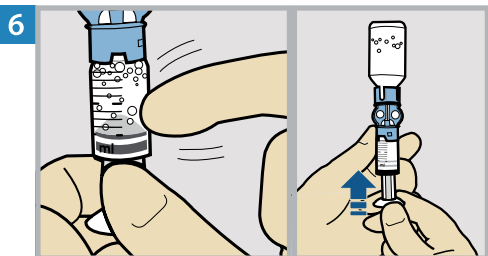


Tap the reservoir to move air bubbles to the top of the reservoir. Push the plunger up to move air into the vial.

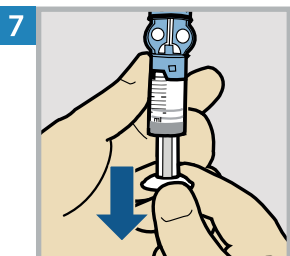

Pull the plunger back down to allow the reservoir to fill with the amount of insulin needed for 2-3 days.


8

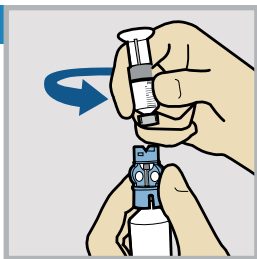

To avoid getting insulin on the top of the reservoir, **flip the vial over again so the reservoir is on top.** Hold the transfer guard and turn the reservoir counter-clockwise and remove the reservoir from the transfer guard.


CONNECT THE INFUSION SET TO THE FILLED RESERVOIR

1 MiniMed™ Quick-set™ Infusion Set

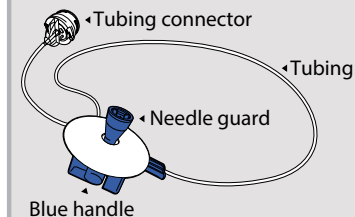


Remove the infusion set tubing from the package. Remove the paper that holds the tubing together and unwind the tubing.

2

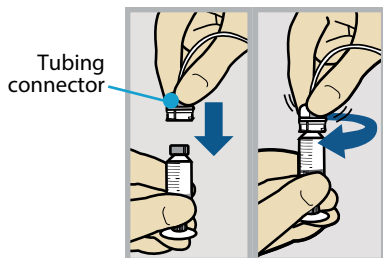


Gently push the tubing connector onto the reservoir. Turn the connector clockwise until it is locked into place.

3

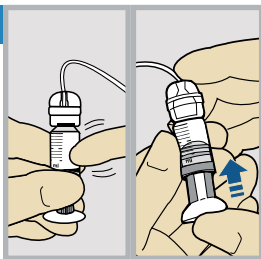

Tap the reservoir to move any air bubbles to the top. Push the plunger slightly to move the bubbles into the tubing.

4

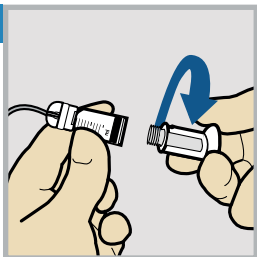

Twist the plunger counter-clockwise to loosen it and to remove it.

PLACE THE RESERVOIR INTO THE PUMP

Note: If the backlight turns off, press any button to turn the screen back on. If you see the Home screen with a red banner, press  to go to the menu screen, and then select . Go through the steps until you reach the screen shown below.

1 New Reservoir & Set

1. Fill new reservoir.
 2. Connect new tubing to reservoir.
- DO NOT CONNECT TO BODY.


Next

Select **Next**.

2 New Reservoir & Set

3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

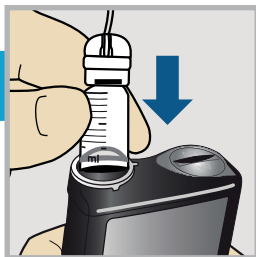


Place the reservoir into the pump.

Do not connect the infusion set to the body.

WARNING: Always confirm that the infusion set is disconnected from the body before placing the reservoir into the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

3

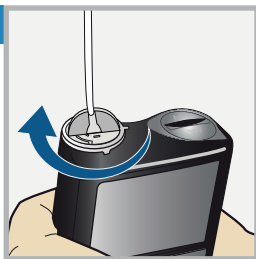

New Reservoir & Set

3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

Turn the reservoir clockwise until the reservoir locks into place, and select **Next**.

LOAD THE RESERVOIR AND FILL THE TUBING

1 Load Reservoir

Select Load and hold until complete.
DO NOT
CONNECT TO
BODY.



Load

Next

Select **Load** and hold  until the checkmark appears on the screen. **Do not** connect the infusion set to the body.

2 Load Reservoir

Complete

DO NOT
CONNECT TO
BODY.



Load

Next

When the checkmark appears, select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before loading the reservoir and filling the tubing. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.


3 Fill Tubing

DO NOT CONNECT TO
BODY.
Hold Fill until drops appear.
Then select Next.

0.0_U

Fill

Next

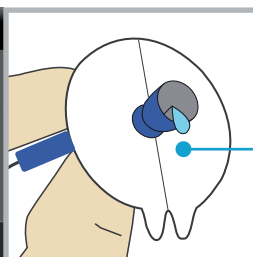

Drops at end
of tubing

Select **Fill** and keep holding  until there are no air bubbles visible in the tubing, and there are drops at the end of the tubing.

Do not connect the infusion set to the body.

4 Fill Tubing

DO NOT CONNECT TO
BODY.


Hold Fill until drops appear.
Then select Next.

11.3_U

Fill

Next



After drops appear, press  and select **Next**.

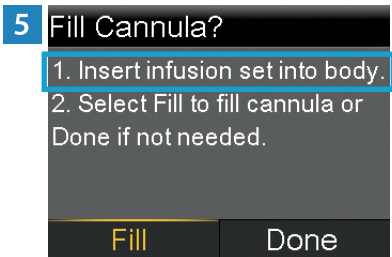


Follow the next steps to insert the infusion set before filling the cannula.

INSERT THE INFUSION SET


1 MiniMed™ Quick-serter™ Insertion Device

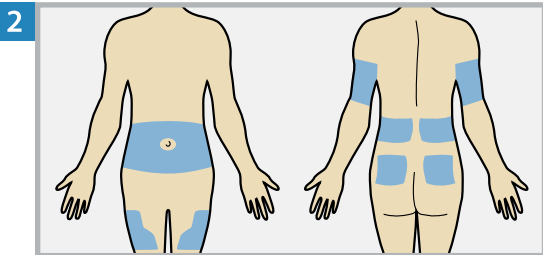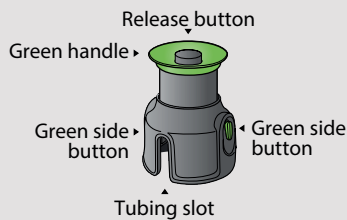

Choose an insertion site from the shaded areas. Clean the insertion site with alcohol or other antiseptic as directed by a healthcare professional.

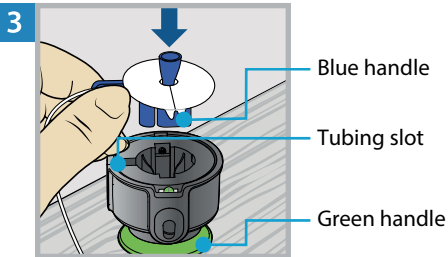

Place the MiniMed Quick-serter insertion device onto a sturdy, flat surface with the green handle facing down. Line up the blue handle of the infusion set with the tubing slot and place the blue handle in the Quick-serter.

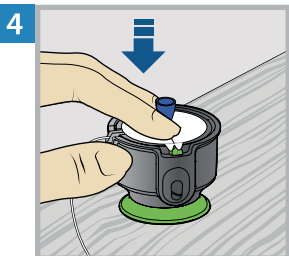


Use two fingers to seat the infusion set inside theserter securely and gently push the set down.




Needle guard

Holding the needle guard, peel the paper from the adhesive on both sides of the needle guard.

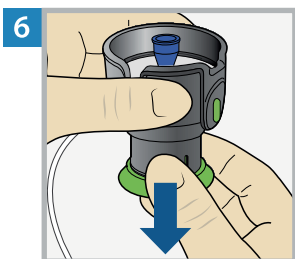


Pull the green handle down until there is a click.

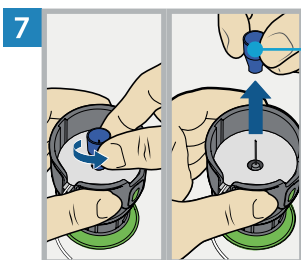


Needle guard

Turn the needle guard to loosen it and pull to remove it.

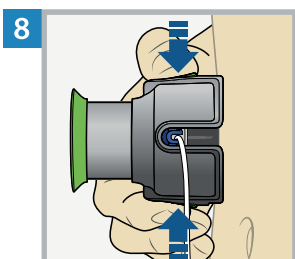


Hold the serter against the cleaned insertion site. Place fingers on the two green side buttons, and press the buttons at the same time. Keep the serter against the cleaned insertion site.

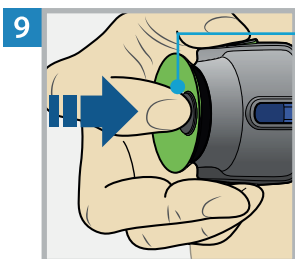


Release button


Press down the release button to release the infusion set from the serter.

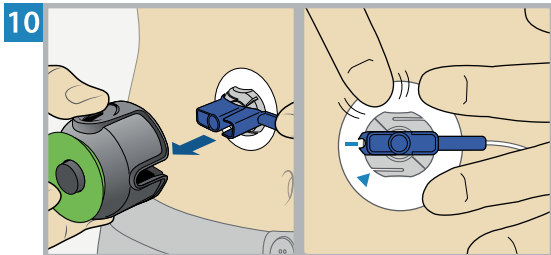

Pull the inserter away from the body. Press the adhesive against the skin all the way around the site.

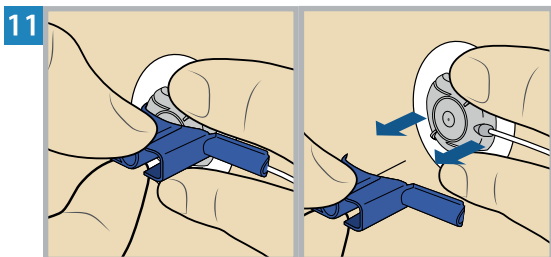


Holding the infusion set, pull the blue handle straight out to remove the needle.

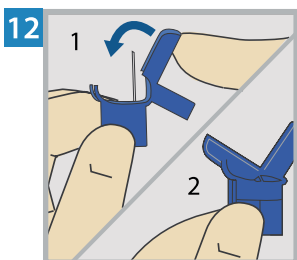

Fold the blue handle over the needle until it locks.

FILL THE CANNULA

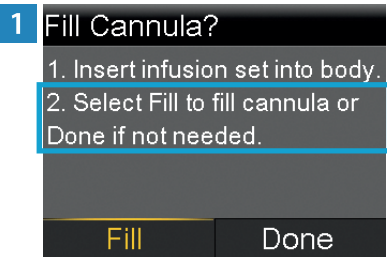

After the infusion set is inserted into the body, select **Fill**.


Note: Always verify that the amount shown in the **Fill amount** field is correct. The pump will remember the fill amount last used. Change the **Fill amount** if needed.



Select **Fill amount** and enter:


- 0.300 U if using the 6 mm cannula
- 0.500 U if using the 9 mm cannula

After entering the cannula fill amount, press .

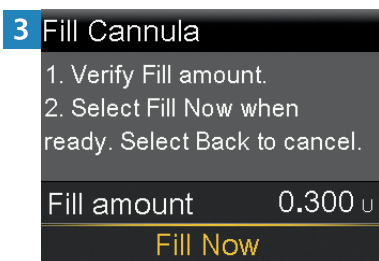

Select **Fill Now**.

The Home screen displays the insulin amount as insulin fills the cannula.

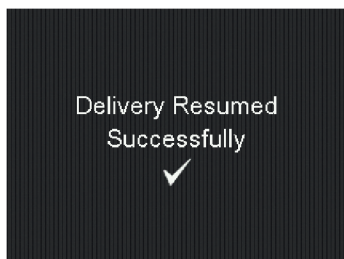

The reservoir and infusion set change is now complete.

Always check blood glucose using a blood glucose meter one to three hours after changing the infusion set or reservoir.


Note: Remove the transfer guard from the vial. Dispose of the transfer guard, the used reservoir, the blue handle and all parts of the used infusion set in an appropriate sharps container and in accordance with local laws.

Disconnecting and Reconnecting the MiniMed Quick-set Infusion Set


DISCONNECT THE INFUSION SET

Follow these steps to disconnect the infusion set. Do not remove the adhesive with the cannula housing from the body.


1 Press the sides of the site connector.

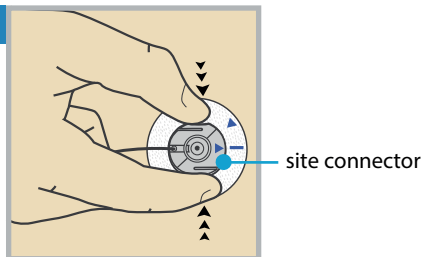

2 Turn the site connector counterclockwise while pressing the sides until the two arrows are aligned.

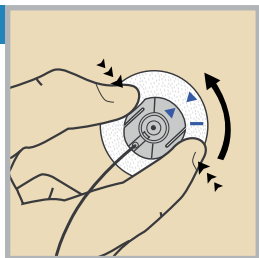


3 Lift up the site connector to remove it from the cannula housing.

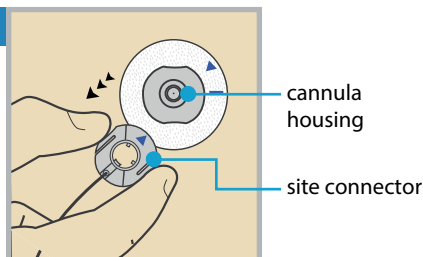

RECONNECT THE INFUSION SET


Follow these steps to reconnect the infusion set.

1

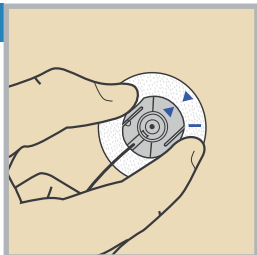

Place the site connector over the cannula housing. Align the blue arrows.

2

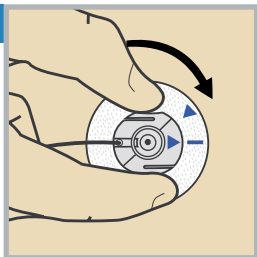

Turn the site connector clockwise without pressing the sides until it locks. The site connector is locked when the arrow on the connector lines up with the blue line on the adhesive.

Changing the MiniMed Reservoir Only



The MiniMed reservoir can be used for up to three days. If the reservoir runs out of insulin and the infusion set has not been used for three full days, the New Reservoir Only option may be used to change the reservoir. The infusion set must still be replaced within three days of use.

Note: The MiniMed Quick-set infusion set and reservoir can be used for up to three days. To avoid wasting insulin, fill the reservoir with enough insulin to last until the next infusion set change.

START HERE

Follow these steps to change the reservoir and keep the current infusion set.

1



Wash hands with soap and water. On the pump, press  to go to the Menu screen.

2

Reservoir & Set
New Reservoir & Set
New Reservoir Only
New Set Only
Fill Cannula


Select , and then select **New Reservoir Only**.

3

New Reservoir Only

Disconnect tubing
from body.

Next

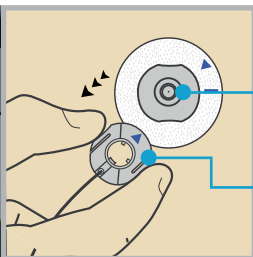

cannula
housing

site connector

Disconnect the tubing from the cannula housing. Do not remove the cannula housing from the body. Select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before rewinding the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.


4 New Reservoir Only

Remove reservoir from pump.
Then select Rewind.

Rewind



Remove the used reservoir from the pump.

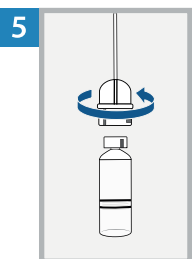

Disconnect the infusion set tubing connector from the used reservoir. Place the tubing on a clean surface. Do not throw the tubing away.

6 New Reservoir Only

Remove reservoir from pump.
Then select Rewind.

Rewind

Select **Rewind**.
Do not connect the infusion set to the body.




7 New Reservoir Only

1. Fill new reservoir.
2. Connect existing tubing to reservoir.
DO NOT CONNECT TO BODY.


Next

Follow the next steps to fill the new reservoir with insulin and to connect the infusion set tubing to the reservoir.
Do not select **Next**.


FILL THE RESERVOIR

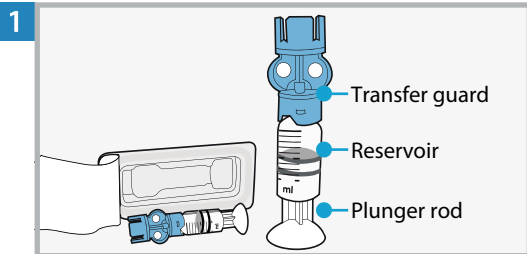


Remove the reservoir from the package. Make sure the insulin vial is at room temperature to reduce the risk of air bubbles.

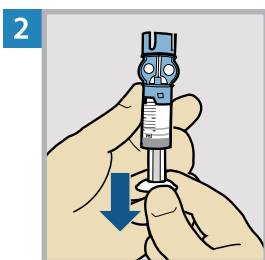

Pull the plunger down based on the planned 2-3 day insulin fill amount.


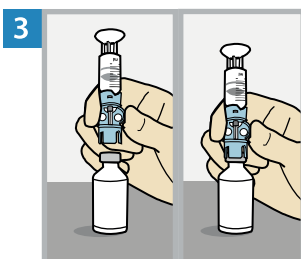


Wipe the top of the vial with alcohol. Place the vial on a sturdy, flat surface. Firmly press the transfer guard onto the vial.

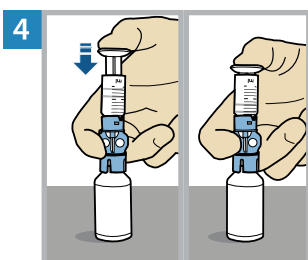


Push and hold the plunger down.


WARNING: Do not use the reservoir or infusion set if insulin or any liquid gets on the top of the reservoir or inside the tubing connector, as shown in the image. Insulin or any liquid can temporarily block the vents. This may result in the delivery of too little or too much insulin, which can cause hyperglycemia or hypoglycemia. If insulin or any liquid gets on the top of the reservoir or inside the tubing connector, start over with a new reservoir and infusion set.

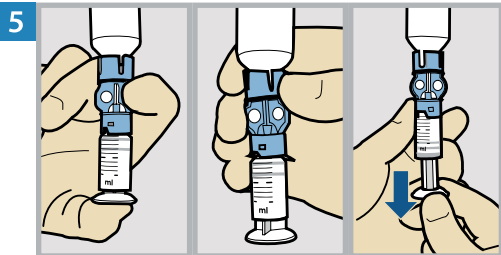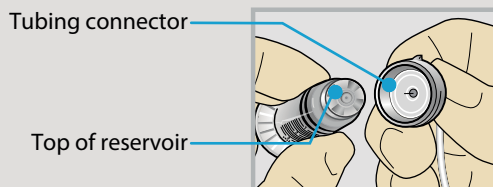

Keeping a thumb on the plunger, flip the vial over so the vial is on top. Release the thumb and pull the plunger down to fill the reservoir with insulin.

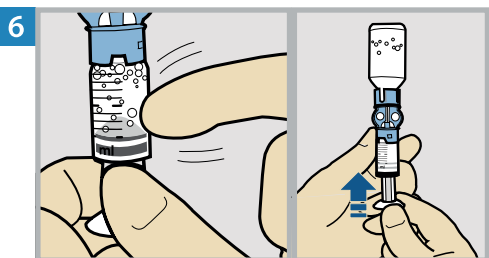


Tap the reservoir to move air bubbles to the top of the reservoir. Push the plunger up to move air into the vial.

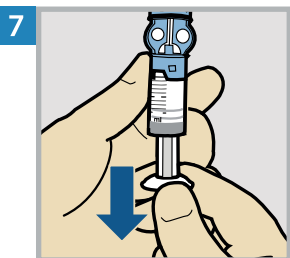

Pull the plunger back down to allow the reservoir to fill with the amount of insulin needed for 2-3 days.

8

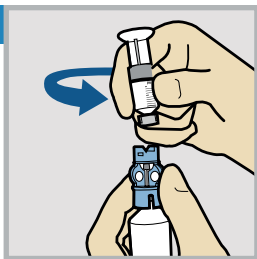


To avoid getting insulin on the top of the reservoir, **flip the vial over again so the reservoir is on top.** Hold the transfer guard and turn the reservoir counter-clockwise and remove the reservoir from the transfer guard.

CONNECT THE INFUSION SET TO THE FILLED RESERVOIR


1

Tubing connector

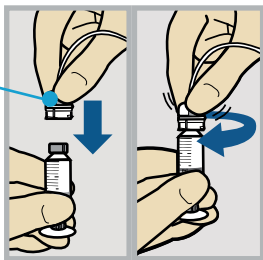


Gently push the tubing connector onto the reservoir. Turn the connector clockwise until it is locked into place.

2

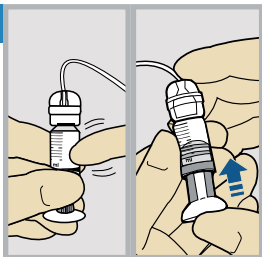

Tap the reservoir to move any air bubbles to the top. Push the plunger slightly to move the bubbles into the tubing.

3

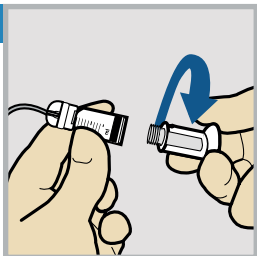

Twist the plunger counter-clockwise to loosen it and to remove it.

PLACE THE RESERVOIR INTO THE PUMP

1 New Reservoir Only

1. Fill new reservoir.
 2. Connect existing tubing to reservoir.
- DO NOT CONNECT TO BODY.


Next

Select **Next**.

2 New Reservoir Only

3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

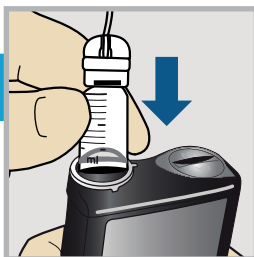


Place the reservoir into the pump.

Do not connect the infusion set to the body.

WARNING: Always confirm that the infusion set is disconnected from the body before placing the reservoir into the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

3

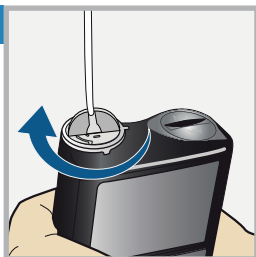

New Reservoir Only

3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

Turn the reservoir clockwise until the reservoir locks into place, and select **Next**.

LOAD THE RESERVOIR AND FILL THE TUBING

1 Load Reservoir

Select Load and hold until complete.
DO NOT
CONNECT TO
BODY.



Load

Next

Select **Load** and hold  until the checkmark appears on the screen. **Do not** connect the infusion set to the body.

2 Load Reservoir

Complete

DO NOT
CONNECT TO
BODY.



Load

Next

When the checkmark appears, select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before loading the reservoir and filling the tubing. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.


3 Fill Tubing

DO NOT CONNECT TO
BODY.
Hold Fill until drops appear.
Then select Next.

0.0_U

Fill

Next

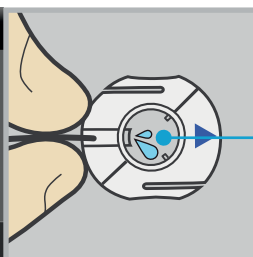

Drops at end
of tubing

Select **Fill** and keep holding  until there are no air bubbles visible in the tubing, and there are drops at the end of the tubing.

Do not connect the infusion set to the body.

4 Fill Tubing

DO NOT CONNECT TO
BODY.

Hold Fill until drops appear.
Then select Next.

11.3_U

Fill

Next


After drops appear, press  and select **Next**.

RECONNECT THE TUBING TO THE INFUSION SITE

1 Reconnect Tubing

1. Connect tubing to the infusion site.
2. Select Resume Basal.

Resume Basal

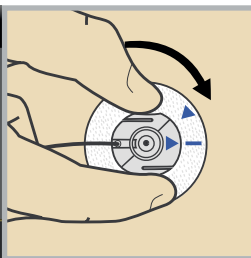


Place the site connector over the cannula housing. Align the blue arrows. Turn the site connector clockwise without pressing the sides until it locks. The site connector is locked when the arrow on the connector lines up with the blue line on the adhesive.

2 Reconnect Tubing

1. Connect tubing to the infusion site.
2. Select Resume Basal.

Resume Basal

It is not necessary to fill the cannula because the infusion set was not changed. Select **Resume Basal**.

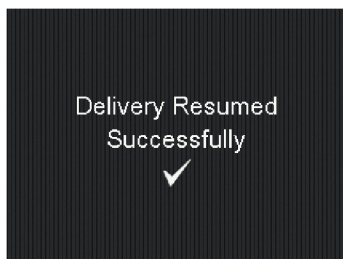


The reservoir change is now complete.

Always check blood glucose using a blood glucose meter one to three hours after changing the infusion set or reservoir.

Note: Remove the transfer guard from the vial. Dispose of the transfer guard and the used reservoir in an appropriate sharps container and in accordance with local laws.

Changing the MiniMed Quick-set Infusion Set Only

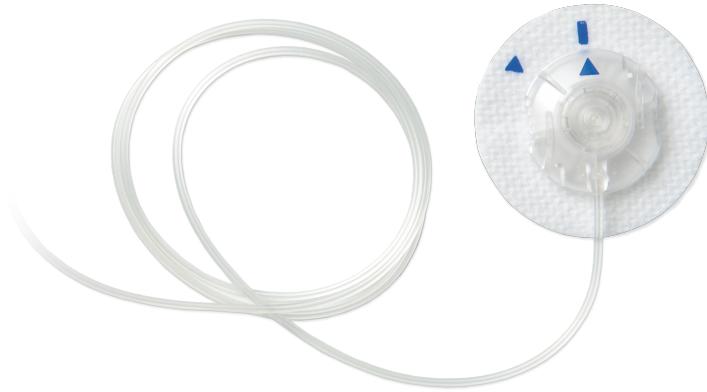

The MiniMed Quick-set infusion set can be used for up to three days. If only the infusion set needs to be changed before three full days of use, the New Set Only option may be used to change the infusion set. The reservoir must still be replaced within three days of use.

START HERE

Follow these steps to change the infusion set and to keep the current reservoir.

1




Wash hands with soap and water. On the pump, press  to go to the Menu screen.

2



Select , and then select **New Set Only**.

3

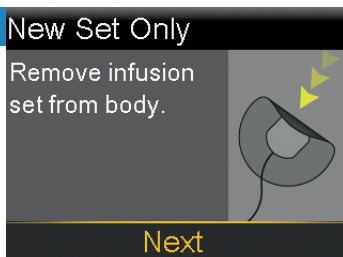


Remove the entire infusion set by loosening the adhesive and pulling the set away from the body. Select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before rewinding the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

4 New Set Only
Remove reservoir from pump.
Then select Rewind.

Rewind

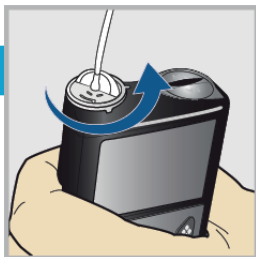

Remove the reservoir from the pump.

5 New Set Only
Remove reservoir from pump.
Then select Rewind.

Rewind

Select **Rewind**.
Do not connect the infusion set to the body.




DISCONNECT THE INFUSION SET FROM THE RESERVOIR

1 New Set Only

1. Disconnect existing tubing from the reservoir.
2. Connect new tubing to existing reservoir.

Next

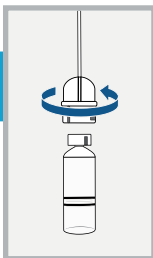

Disconnect the tubing connector of the used infusion set tubing from the reservoir. Place the reservoir on a clean surface. Do not throw the reservoir away.

CONNECT THE INFUSION SET TO THE RESERVOIR


1 New Set Only

1. Disconnect existing tubing from the reservoir.
2. Connect new tubing to existing reservoir.

Next


Follow these steps to connect the infusion set to the reservoir.

2 MiniMed™ Quick-set™ Infusion Set

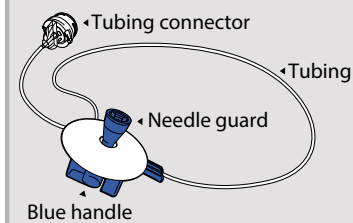


Remove the infusion set tubing from the package. Remove the paper that holds the tubing together and unwind the tubing.

3

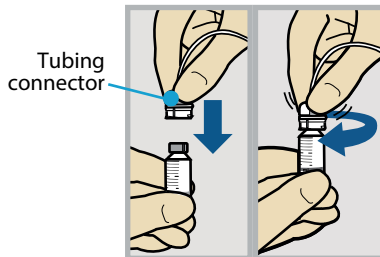

Gently push the tubing connector onto the reservoir. Turn the connector clockwise until it is locked into place.

WARNING: Do not use the reservoir or infusion set if insulin or any liquid gets on the top of the reservoir or inside the tubing connector, as shown in the image. Insulin or any liquid can temporarily block the vents. This may result in the delivery of too little or too much insulin, which can cause hyperglycemia or hypoglycemia. If insulin or any liquid gets on the top of the reservoir or inside the tubing connector, start over with a new reservoir and infusion set.

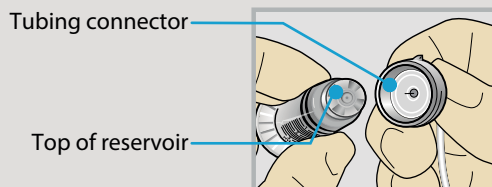

PLACE THE RESERVOIR INTO THE PUMP

1 New Set Only

1. Disconnect existing tubing from the reservoir.
2. Connect new tubing to existing reservoir.


Next

Select **Next**.

2 New Set Only

3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

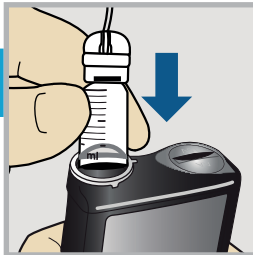


Place the reservoir into the pump.

Do not connect the infusion set to the body.

WARNING: Always confirm that the infusion set is disconnected from the body before placing the reservoir into the pump. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.

3

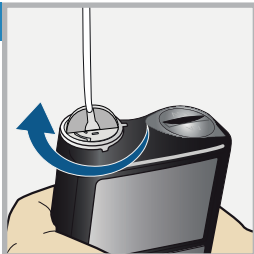

- New Set Only
3. Place reservoir into pump and lock.
- DO NOT CONNECT TO BODY.

Next

Turn the reservoir clockwise until the reservoir locks into place, and select **Next**.

LOAD THE RESERVOIR AND FILL THE TUBING

1 Load Reservoir

Select Load and hold until complete.
DO NOT
CONNECT TO
BODY.



Load

Next

Select **Load** and hold  until the checkmark appears on the screen. **Do not** connect the infusion set to the body.

2 Load Reservoir

Complete

DO NOT
CONNECT TO
BODY.



Load

Next

When the checkmark appears, select **Next**.

WARNING: Always confirm that the infusion set is disconnected from the body before loading the reservoir and filling the tubing. Failing to disconnect the infusion set from the body may result in an accidental infusion of insulin, and may cause hypoglycemia.


3 Fill Tubing

DO NOT CONNECT TO
BODY.
Hold Fill until drops appear.
Then select Next.

0.0_U

Fill

Next

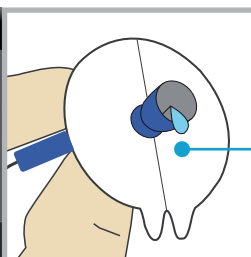

Drops at end
of tubing

Select **Fill** and keep holding  until there are no air bubbles visible in the tubing, and there are drops at the end of the tubing.

Do not connect the infusion set to the body.

4 Fill Tubing

DO NOT CONNECT TO
BODY.

Hold Fill until drops appear.
Then select Next.

11.3_U

Fill

Next

After drops appear, press  and select **Next**.

5 Fill Cannula?

1. Insert infusion set into body.
2. Select Fill to fill cannula or Done if not needed.



Fill

Done


Follow the next steps to insert the infusion set before filling the cannula.

INSERT THE INFUSION SET


1 MiniMed™ Quick-serter™ Insertion Device

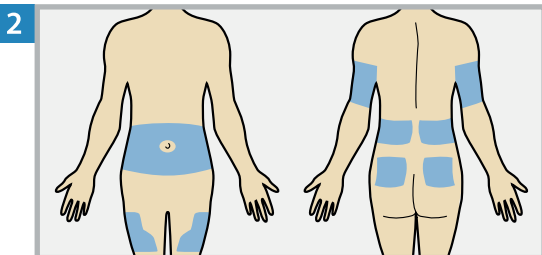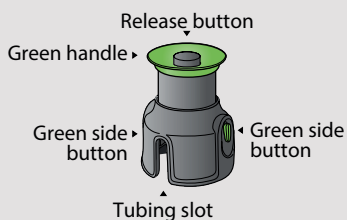

Choose an insertion site from the shaded areas. Clean the insertion site with alcohol or other antiseptic as directed by a healthcare professional.

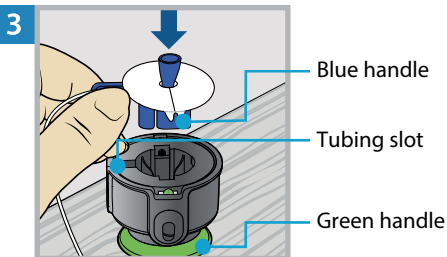

Place the MiniMed Quick-serter insertion device onto a sturdy, flat surface with the green handle facing down. Line up the blue handle of the infusion set with the tubing slot and place the blue handle in the Quick-serter.

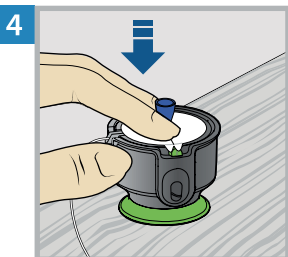

Use two fingers to seat the infusion set inside the serter securely and gently push the set down.


5



Needle guard


Holding the needle guard, peel the paper from the adhesive on both sides of the needle guard.

6

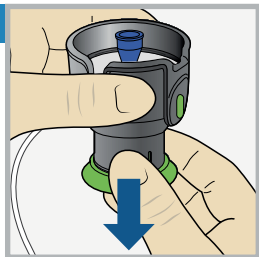

Pull the green handle down until there is a click.


7

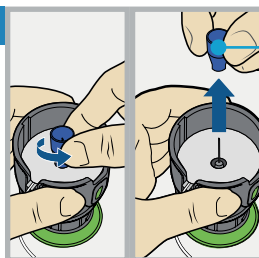

Needle guard


Turn the needle guard to loosen it and pull to remove it.

8

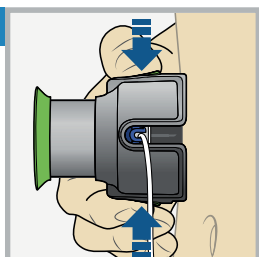


Hold the inserter against the cleaned insertion site. Place fingers on the two green side buttons, and press the buttons at the same time. Keep the inserter against the cleaned insertion site.

9

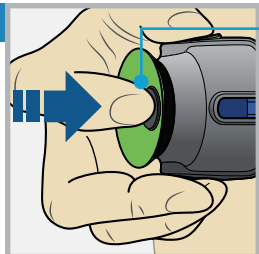


Release button


Press down the release button to release the infusion set from the inserter.

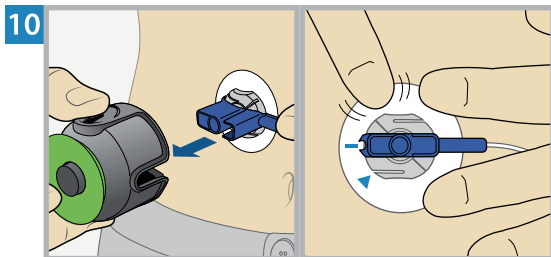

Pull the inserter away from the body. Press the adhesive against the skin all the way around the site.

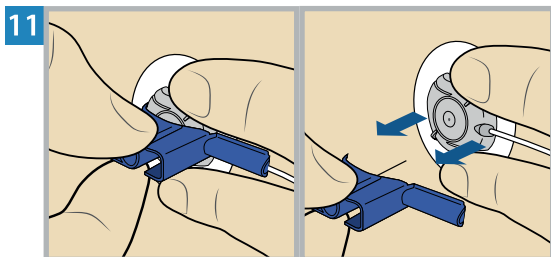

Holding the infusion set, pull the blue handle straight out to remove the needle.

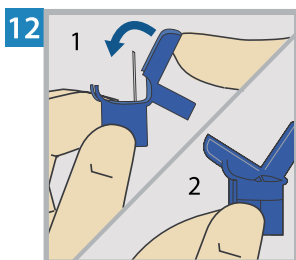

Fold the blue handle over the needle until it locks.

FILL THE CANNULA

1 Fill Cannula?

1. Insert infusion set into body.
2. Select Fill to fill cannula or Done if not needed.

Fill

Done

After the infusion set is inserted into the body, select **Fill**.

Note: Always verify that the amount shown in the **Fill amount** field is correct. The pump will remember the fill amount last used. Change the **Fill amount** if needed.

2 Fill Cannula

1. Verify Fill amount.
2. Select Fill Now when ready. Select Back to cancel.

Fill amount

--- U

Fill Now

Select **Fill amount** and enter:

- 0.300 U if using the 6 mm cannula
- 0.500 U if using the 9 mm cannula

After entering the cannula fill amount, press .

3 Fill Cannula

1. Verify Fill amount.
2. Select Fill Now when ready. Select Back to cancel.

Fill amount

0.300 U

Fill Now

Select **Fill Now**.

The Home screen displays the insulin amount as insulin fills the cannula.

Delivery Resumed
Successfully



The infusion set change is now complete.

Always check blood glucose using a blood glucose meter one to three hours after changing the infusion set or reservoir.

Note: Dispose of the blue handle and all parts of the used infusion set in an appropriate sharps container and in accordance with local laws.

Inserting the Guardian™ 4 Sensor with One-press Serter




INSERT THE SENSOR


Choose an insertion site that has an adequate amount of subcutaneous fat. The Guardian 4 sensor has been studied and is approved for use in the following sensor insertion sites by persons of the following ages:

Approved Age	Sensor Insertion Site
7-17	Arm
18 and older	Abdomen and Arm

Ages 7-17 years

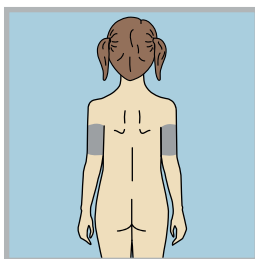


Back of the upper arm

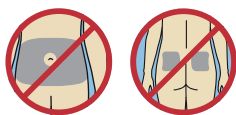


Insertion on abdomen and upper buttocks for ages 7-17 years is not approved.

Note: Assistance from another person may be needed for sensor insertion into the back of the upper arm. If assistance is not needed, a mirror may be helpful for self-insertion.

Ages 18 years and older

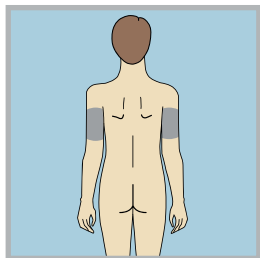

**Back of the
upper arm**

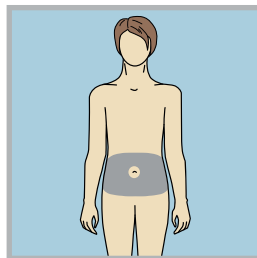

Abdomen



Insertion on upper buttocks for ages 18 years and older is not approved.

INSERT A NEW SENSOR


Wash hands with soap and water and clean the insertion site with alcohol.

1




Open the sensor package. Pull the corner of the paper covering to open the sensor package.

2a

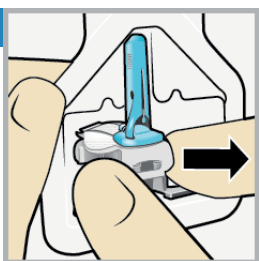

Hold the sensor by the plastic pedestal. Remove the sensor with attached pedestal from the packaging by holding the pedestal only. Place the sensor with the pedestal on a clean, flat surface (such as a table).


2b

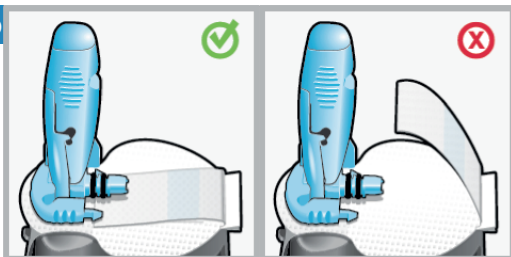

Correct

Incorrect

Tuck the adhesive tab. Make sure that the sensor's adhesive tab is tucked under the sensor connector and snaps.

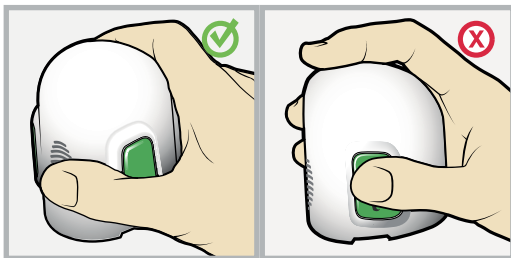


Correct

Incorrect


Note: Refer to the illustrations for correct and incorrect ways to hold the serter for loading.

3

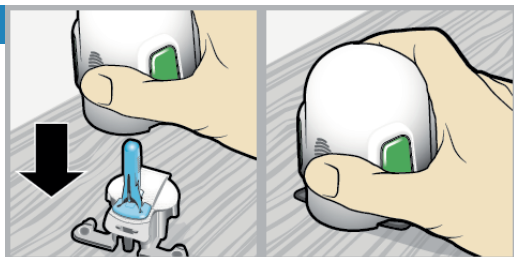

Load the sensor into the serter. Grip the serter exactly as shown **with a thumb on the thumb print on the serter. Do not hold the side buttons.** Push the serter down onto the pedestal until the base of the serter sits flat on the table.

4

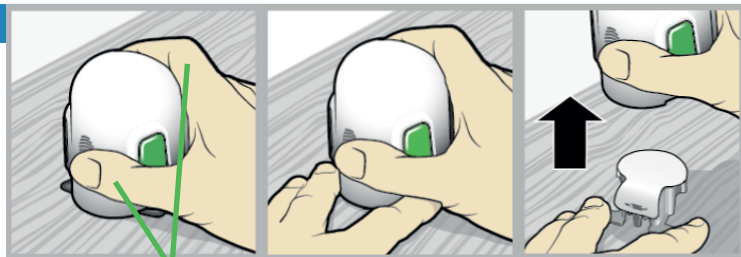

Fingers are **NOT** holding the side buttons.

Detach the serter from the pedestal. To detach the serter from the pedestal, grip the serter as shown with a thumb on the thumb print on the serter. With the other hand, place two fingers on the pedestal arms, and slowly pull the serter straight up.

Note: Make sure that the pedestal is firmly on the table before pulling the serter away.


CAUTION: Do not detach the pedestal from the serter in mid-air as this may damage the sensor.

5a




Place the serter on the body. Hold the serter steadily against the cleaned insertion site without pushing the serter too deeply into the skin.


Note: Failing to hold the serter securely flat against the body may allow the serter to spring back after pressing buttons and result in improper insertion of the sensor.

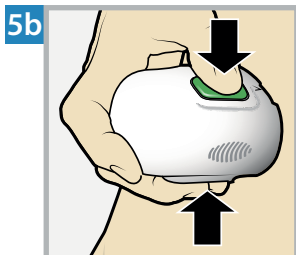


5b Insert the sensor. Press the **bump** on both buttons and release them at the same time.

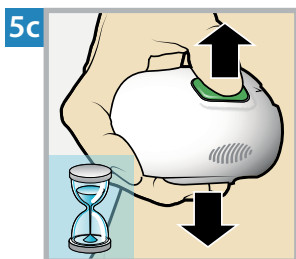

5c Hold the serter against the body. Continue holding the serter against the body for at least five seconds to allow the adhesive to stick to the skin.

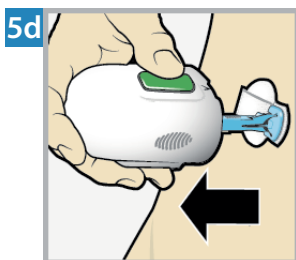


5d Remove the serter from body. Slowly pull the serter away from the skin, making sure the buttons are not pressed.

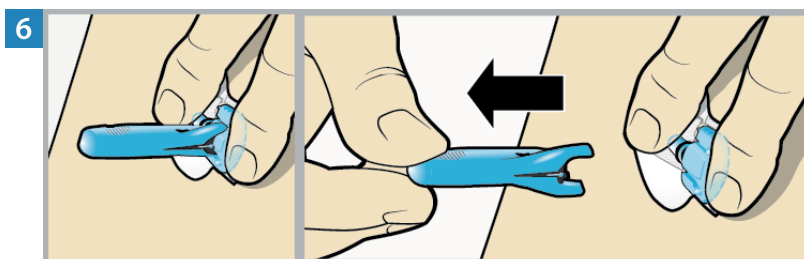


6 Remove the needle housing. Gently hold the base of the sensor against the skin with one hand. With the other hand, hold the needle housing **at the top** and slowly pull it straight away from the sensor. Dispose of the needle housing in a sharps container.

Note: An optional, liquid adhesive, such as Skin Tac™*, may be used on the skin. Apply the adhesive and allow the skin to dry before removing the adhesive pad liner.

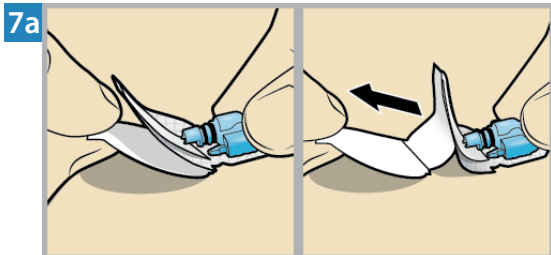


Remove the adhesive pad liner. Hold the sensor in place and gently remove the liner from under the adhesive pad.

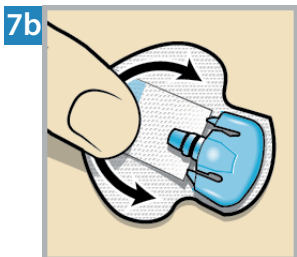


Press the entire adhesive pad to the skin. Firmly press the adhesive pad against the skin and smooth the entire adhesive pad so it sticks to the skin.

Note: The sensor adhesive is sensitive to pressure. Continue applying pressure on the adhesive to ensure that the sensor remains inserted in the skin for 7 days of wear.

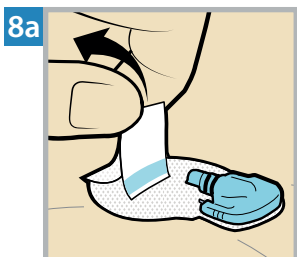


Untuck the adhesive tab. Untuck the adhesive tab from under the sensor connector.

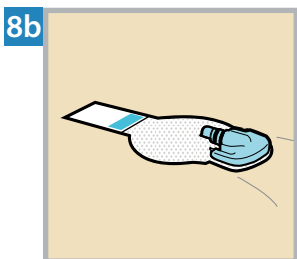

Straighten the adhesive tab. Straighten the adhesive tab so it lies flat against the skin, but do not remove the adhesive liner yet.


TAPE THE SENSOR

1

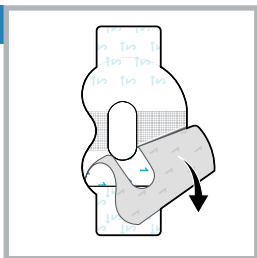

Remove the liner marked 1.

2

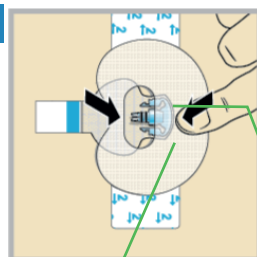


Apply the tape as shown and press down firmly.

Both sensor and skin are taped.

Wide part of tape covers half of sensor base.


3

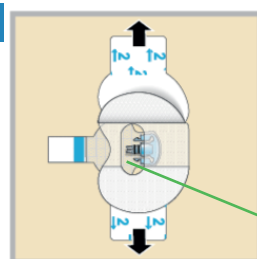

Remove the liner marked 2 from each side.

Connector and snaps in hole of tape.


4

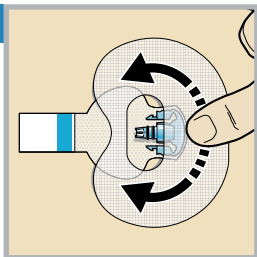

Smooth the tape.

CONNECT THE TRANSMITTER


1

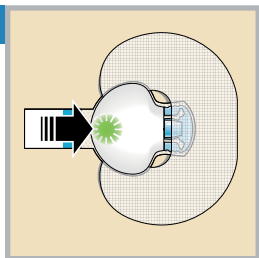

Connect the transmitter to the sensor.

Note: Wait for the green light on the transmitter to flash. If the green light does not flash, refer to Troubleshooting section of the transmitter user guide.


2

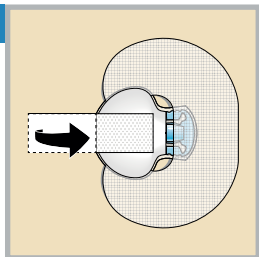


Cover the transmitter with the adhesive tab.

Note: Do not pull the tab too tightly.

3

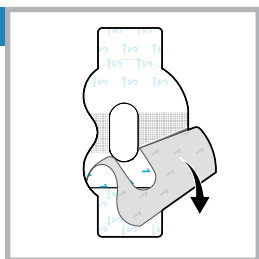


To apply the 2nd tape, remove the liner marked 1.

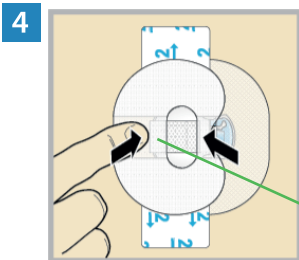


Wide part of tape covers end of transmitter and skin.

Rotate the 2nd tape and place the tape over the transmitter. Press down firmly.

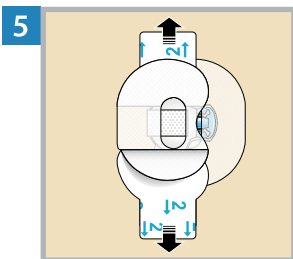


Remove the liner marked 2 from each side.

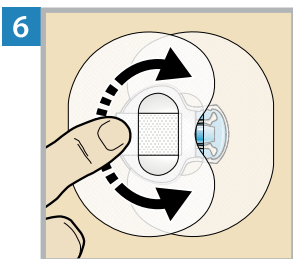

Smooth the tape.

Note: Check the site regularly. Apply additional off-the-shelf tape if the sensor and the transmitter are not secure.

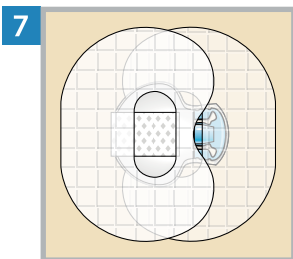

This image is an example of the tape applied correctly.

START THE SENSOR

1 Sensor connected

If new sensor, select Start New. If not, select Reconnect.

After **Sensor connected** message appears, press ▼. The message typically takes less than a minute to appear, but may take up to 5 minutes.

2 Sensor connected

Reconnect.

Start New Sensor

Reconnect Sensor

Select **Start New Sensor**.

3

09:00

Sensor warm up 2:00 hr

0.425 U/hr • 0.1 U Act. Insulin

Sensor warm up X:XX hr will appear on the Home screen. After the warm up is complete, the pump begins receiving SG readings.

Medtronic



Medtronic MiniMed

18000 Devonshire Street

Northridge, CA 91325

USA

1 800 646 4633

+1 818 576 5555

www.medtronicdiabetes.com

Medtronic Canada ULC

99 Hereford Street

Brampton ON L6Y 0R3

800 284 4416

www.medtronicdiabetes.ca



(01)00199150077441(10)NA

© 2026 Medtronic
M080724C001_1
2026-04-09



M080724C001