


Wallet Card

Medtronic

MiniMed™ 780G



Main Menu

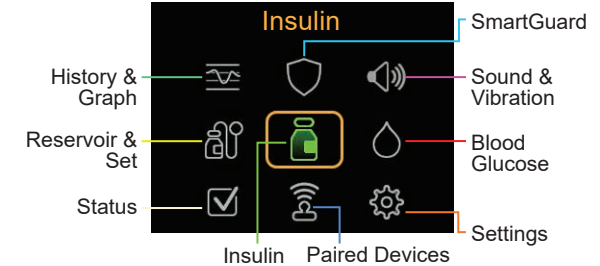

History & Graph History Sensor Glucose Review* Graph Time in Range*	SmartGuard SmartGuard Checklist Temp Target* SmartGuard Settings* SmartGuard*
Sound & Vibration Silence Sensor Alerts* Volume Sound Vibration Alert Settings	Reservoir & Set New Reservoir & Set New Reservoir Only New Set Only Fill Cannula*
Insulin Bolus Basal Suspend All Delivery Delivery Settings	Blood Glucose Enter BG
Status Suspend All Delivery SmartGuard Checklist Pump Sensor*	Paired Devices Pair New Device Pair CareLink Mobile** Meter** CGM**
Settings Alert Settings Delivery Settings Device Settings	<i>*In some cases, menu choices appear, but are greyed out if certain settings are turned on or off.</i> <i>**Menu choices appear when the device is paired to the pump.</i>

status icons

09:00

current time

6.1

mmol/L

BG reading

Active Insulin

1.0 U

active insulin

status icons

7 8 9

20 17 14 11 8 5 2.2

sensor graph

09:00

current time

9.2

mmol/L

SG reading

0.5 U

Act. Insulin


active insulin

suspend by sensor icon

trend arrows

Home screen in SmartGuard™ feature

status icons

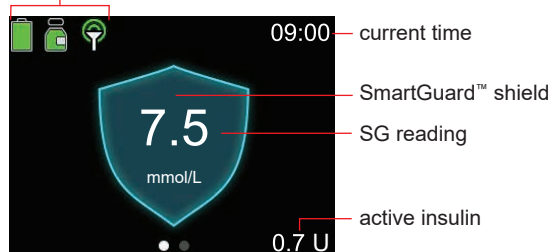

Stopping and resuming insulin delivery

Stopping all insulin delivery

The Suspend All Delivery option stops all insulin delivery (basal and bolus).

To stop all insulin delivery:

1. From the Home screen, press  and select .
2. Select **Suspend All Delivery**.
3. Select **Yes** when prompted.

Resuming basal insulin delivery

The Resume Basal option restarts the basal delivery after the pump has been suspended. The Resume Basal option does not restart a bolus delivery.

To resume basal insulin delivery:

1. From the Home screen, press  and select .
 2. Select **Resume Basal**.
 3. Select **Yes** when prompted.
- NOTE:** Basal insulin delivery can be manually resumed during a Suspend before low or Suspend on low event.

Stopping a bolus delivery

The Stop Bolus option stops a bolus delivery only. The Stop Bolus option does not stop the basal delivery.

To stop a normal bolus, or the Now portion of a Dual Wave™ bolus:

1. From the Home screen, press  and select .
2. Select **Stop Bolus**.
3. Select **Yes** when prompted.
4. Select **Done**.

To stop a Square Wave™ bolus or the Square portion of a Dual Wave™ bolus:

1. From the Home screen, press  and select .
2. Select **Bolus**.
3. Select **Stop Bolus**.
4. Select **Yes** when prompted.
5. Select **Done**.

Manually resuming basal insulin delivery during a Suspend before low or Suspend on Low event.

If unable to wait for the pump to automatically resume the basal insulin during a Suspend before low or Suspend on low event, manually resume the basal delivery.

To manually resume basal insulin delivery:

1. From the Home screen, press  and select .
2. Select **Resume Basal**.
3. Select **Resume Basal** again.
4. Select **Yes** when prompted.

© 2023 Medtronic. Medtronic and Medtronic logo are trademarks of Medtronic.

™* Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

MiniMed™, SmartGuard™, Dual Wave™, and Square Wave™ are trademarks of Medtronic MiniMed, Inc.

Medtronic



Medtronic MiniMed
18000 Devonshire Street
Northridge, CA 91325 USA
1800 646 4633
+1818 576 5555
www.medtronicdiabetes.com

EC REP

Medtronic B.V.
Earl Bakkenstraat 10
6422 PJ Heerlen
The Netherlands

C €0459

2023-11-14

M028312C001_2